

Compat ella®

Compat Ella® Pump Combiset 1L ENFit™

The Compat Ella® Pump Combiset 1L is intended for administration of enteral nutrition, fluids or medication and is suitable for use with the Compat Ella® feeding pump. NON-Sterile.

7eatures at a glance

Built-in Safety: All Compat Ella® administration sets feature patented fail-safe Auto-Anti-Free-Flow (AAFF) technology. AAFF is incorporated into the cassette and is always automatically activated to ensure patient safety.


Safer Connectors: The ENFit[™] connector is designed to be IV INcompatible to promote patient safety by helping to avoid accidental misconnections.

Convenience: The integrated 1L reservoir offers a convenient open system with a built-in hanger and screw cap. The large top-fill opening allows to refill the bag with nutrition or water without stopping the feeding program, or removing the reservoir from the hanger.

Food Safe Materials: All Compat Ella® administration sets are DEHP, BPA and Latex-Free*. The Flexibaggle® 1L reservoir is made of Ethylene-vinyl Acetate (EVA), a soft and food friendly material.

Convenient:


- Enteral feeding tubes: ENFit[™], ENLock, Funnel**
- · Medication dispensers: ENFit™, ENLock, Oral, Luer

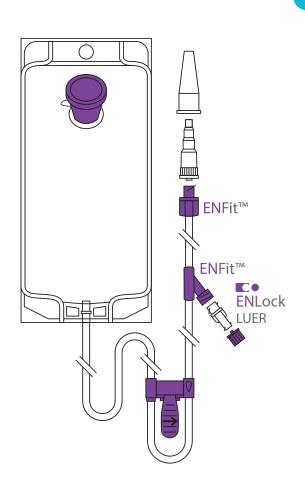

Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella® Pump Combiset 1L ENFit™	12277771	30

Compat Ella® Pump Combiset 1L ENFit™

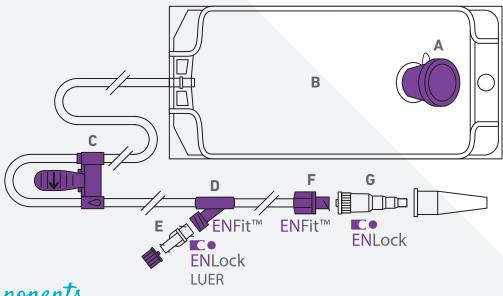

Components

- A Large top-fill opening with built-in hanger and screw cap
- B Flexibaggle® 1L container made of EVA with graduations (at 50 mL and every 100 mL) and tags for patient name and other important information (formula, volume, date and rate)
- C Compat Ella® pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- **D** ENFit[™] Y medication port
- **E** Removable medication port adapter (Luer, Oral and ENLock compatible) and cap
- **F** ENFit[™] distal end (patient side tube connection)
- **G** Removable adapter with cap (ENLock and legacy funnel tubes compatible**)

Specifications

Sterility: NON-sterile.

Packing: Individually packed in cases of 30 pieces (58.5 x 38.5 x 16 cm).

Storage: Store at room temperature. Keep away from sunlight and humidity.

Shelf life: 3 years from production date. Expiry date on the package.

Duration of use: 24 hours. Single Use.

Tube Length: Distance from bag to patient side tube connection: 230 cm.

Distance from Flexibaggle® to pump: 50 cm.

Tube diameters: Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.

CE Mark: Complies with the European Medical Device Directive 93/42 EEC. A CE Certification has been

issued by an accredited Notified Body: TÜV SÜD (CE0123).

Ecological Footprint: The 'Grüne Punkt' symbol signals that Nestlé Health Science assumes its producer responsibi-

lity and contributes to the financing of the recycling.

Materials: Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free*. Additional details on request.

To obtain further information please contact your Nestlé Health Science representative or go to www.compatella.com

*The parts do not exceed the maximum concentrations of 0.1% total DEHP (Di (2-Ethyl Hexyl) Phthalate) content.

No Latex was intentionally added during the manufacturing process.

**Some incompatibilities may exist.

Unless stated otherwise, all trademarks are owned by Société des Produits Nestlé S.A or used with permission. www.compatella.com/patent.
 Legal Manufacturer: Cedic Srl.
 For more information on ENFit™ visit www.stayconnected.com

NHSMD00209 Rev.0 16Feb2016

